

PHILIPS

Home Nebulizer

Compressor
nebulizer system

Affordable high quality aerosol therapy

Building on more than 30 years of experience in aerosol therapy, Philips Home Nebulizer offers high quality aerosol therapy for the effective management of respiratory conditions, at an affordable price.

Reliable and easy to use, Philips Home Nebulizer is suitable for nebulizing most aerosol therapy drugs.



Economical choice

Affordable aerosol therapy - suitable for children and adults...



Effective therapy

Reliable aerosol therapy delivered using the clinically proven SideStream nebulizer technology**.



Treatment in 6 mins*

A compact design that provides delivery of consistent aerosol medication treatments in as little as 6 minutes.



Easy to use

The smallest of the Philips home compressor systems, Philips Home Nebulizer is also lightweight and simple to operate.




* 2.5ml Salbutamol solution drug

** An evaluation of two aerosol delivery systems for rhDNase
- PL Shah et al. European Respiratory Journal 1997

Reliable aerosol delivery you can trust

Philips Home Nebulizer is combined with our highly efficient SideStream nebulizer. SideStream has been designed to draw more air in, via an active venturi system, resulting in fast drug delivery and short treatment times. It also has a unique diamond jet design which helps to minimize wear from frequent use and aids consistent delivery of a high-quality aerosol.

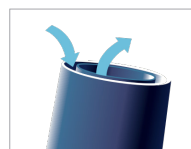

By utilizing a venturi, SideStream produces a flow rate of 16 l/min



Innovative jet design producing consistently high quality aerosol



Product specification table

Dimensions	6.7x9.3x3.8 inches
Weight	2.6lbs (1.17 Kg)
Noise level	60 +/- 4dBA
Mode of operation	Intermittent (30 min ON/30 min OFF)

External packaging – weights and dimensions

Product description	Size	Weight
Philips Home Nebulizer (unit carton) 1 EA	8.8x5.8x6.5 inches	3.32lbs (1.51 Kg)
Philips Home Nebulizer shipped in cases of 6	18.5x9.4x14.2 inches	21.1 lbs (9.6 Kg)

Ordering information 230V/50Hz EU plug

Product description	Reorder number	UOM
Philips Home Nebulizer with SideStream disposable kit, adult and pediatric masks	1130529	1 EA
SideStream disposable kit with adult and pediatric masks	2350	1 case of 40
Replacement filters	1083128	1 pack of 5

Ordering information 120V/60Hz US plug

Product description	Reorder number	UOM	NDC number	HPCPS code
Philips Home Nebulizer with SideStream reusable and disposable nebulizers and mouthpiece	1130530	1 EA	08373-0530-00	E0570
Philips Home Nebulizer with SideStream disposable nebulizer and mouthpiece	1130531	1 EA	08373-0531-00	E0570
SideStream reusable kit	HS860	1 case of 10	N/A	A7005
SideStream disposable kit	HS800	1 case of 50	N/A	A7003
Replacement filters	1083128	1 pack of 5	08373-0528-00	A7014
Optional accessories:				
Adult masks	1025531	1 case of 50	N/A	A7015
Pediatric masks	1025529	1 case of 50	N/A	A7015

Find out more – www.philips.com/respironics

